

NEW – NEW – NEW THE BITE SPLINT WITH THERMO-EFFECT



excited
patients

pleased
dentists

relaxed
technicians

Nothing presses, nothing strains

**the milled PREMIOtemp CLEAR FLEX bite
splint is just super comfortable to wear.**


PREMIOtemp CLEAR FLEX

the **NEW** standard for milled **bite splints!**

- **PREMIOtemp CLEAR FLEX is the new** and very innovative milling blank for **CAD/CAM milled nightguards** to our primotec digital products portfolio.
- The CLEAR FLEX material is a **transparent PMMA acrylic with thermo-effect**. This effect procures that the milled bite splint, when warmed up, automatically adjusts to the dental arch in the patients mouth, without losing the materials surface hardness.
- The thermo-effect is achieved, when the milled PREMIOtemp CLEAR FLEX nightguard is **pre-heated in 40°C to 50°C warm water** prior to being inserted. This way, during insertion, the bite splint perfectly adapts itself in the patient's mouth.

>>>>

CLEAR FLEX

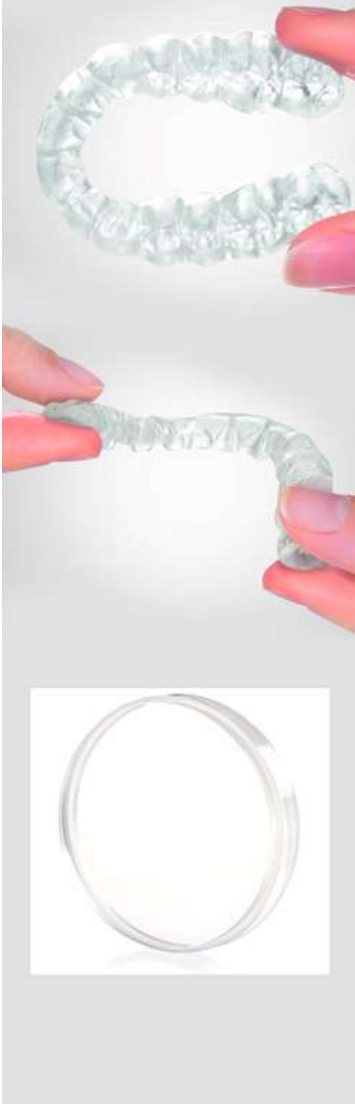

THE NEW STANDARD FOR MILLED BITE SPLINTS



O
N
E
F
R

The result is a milled nightguard that **perfectly fits the dental arch without any tension**. Due to the surface hardness of the material, occlusal contacts and canine- as well as anterior guidance remain where they had been set in the lab.

Especially comfortable for the dentist and the lab is the fact, that PREMIOtemp CLEAR FLEX, due to the described thermo-effect, **can compensate for moderate inaccuracies of the impression**. This way, the material is especially appropriate **for intra oral scanning of the complete dental arch**.

CONCLUSION:

- **excited patients**
- **pleased dentists**
- **relaxed technicians**

due to milled PREMIOtemp CLEAR FLEX nightguards from primotec

Contact: 203-682-6429

or: info@primotecusa.com