

QUICK MODELING DIRECTLY FROM THE SYRINGE

metablue

light cured Metacon system modeling composite

metablue is a light cured modeling composite gel from the syringe,
that burns out cleanly and completely and is aligned specifically
to be used with the metacon light cured wax system.



QUICK MODELING DIRECTLY FROM THE SYRINGE

uses

- individual implant abutments
- secondary parts of combined cases, such as telescopic or conus crowns, interlocks, RSS, millings, individual attachments, etc.
- tension free connection of Metacon implant bars, bridges or partial segments
- smoothing and/or sealing of ground Metacon surfaces
- fixing of retention beads on secondary telescopic or conus crowns
- easy attaching of the Metacon sprues to the modelled case





www.primotecusa.com

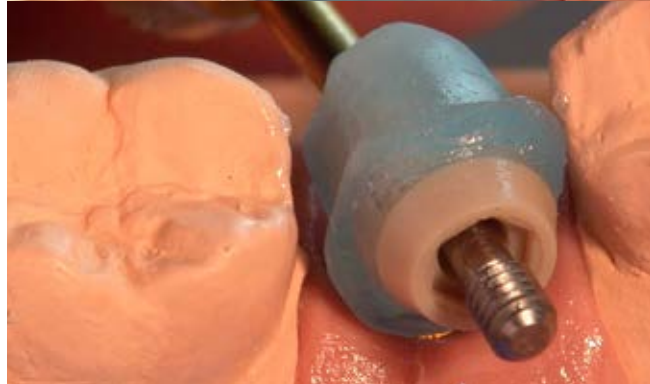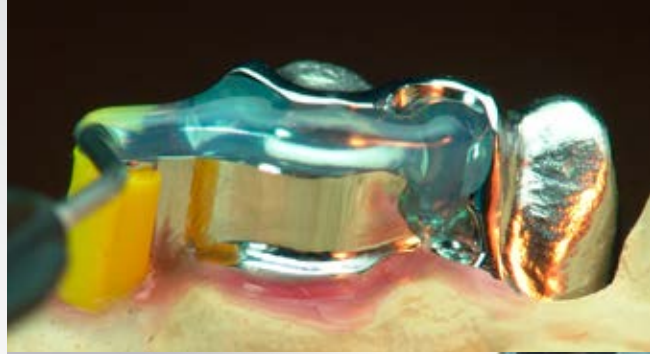


advantages

- ready to use, one component material from the syringe
- burns out cleanly and completely
- very good and precise modeling characteristics
- connects directly to the light cured Metacon wax
- is processed "cold" – no danger of thermal tension when connecting Metacon segments
- no separator needed on zirconium, metal or porcelain surfaces
- no danger of "melting through" on the block-out or relief wax (in case of partials)
- long working time (> 20 min.), short light curing time (< 5 min.)
- absolutely precise after polymerisation – no noticeable shrinkage or distortion
- secondary parts can be lifted off without any problem
- easy to grind with conventional carbide burs or universal rubber polishers

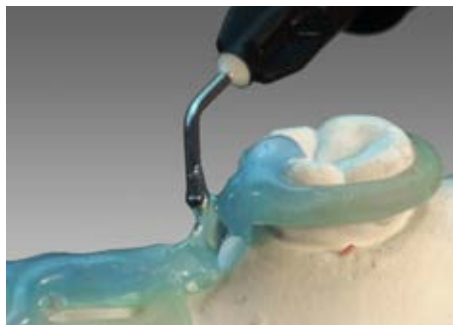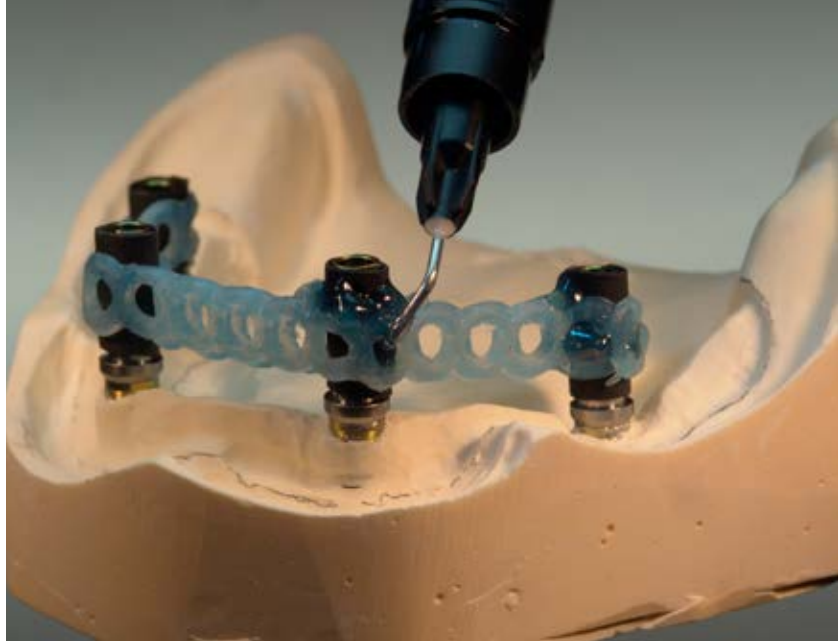

ORDERING INFO

Item #	Item	Price
MB100	 <p>metablue light cured Metacon System modeling composite, 4 syringes at 3 g / each, 12 hollow needle tips</p>	\$ 79.00

quick modeling directly from the syringe

with metablue the metacon system becomes even easier in handling with even more perfect results



Primotec USA
215 Westport Ave.
Norwalk, CT 06851
203-682-6429

info@primotecusa.com
www.primotecusa.com